

Our Most Versatile Laboratory Clamping System



TALBOYS


201 Wolf Drive • PO Box 87 • Thorofare, NJ 08086-0087

Phone: 1-888-TALBOYS • Fax: 856-686-1601

www.talboys.com • E-mail: info@talboys.com

TALBOYS

Ultra Flex

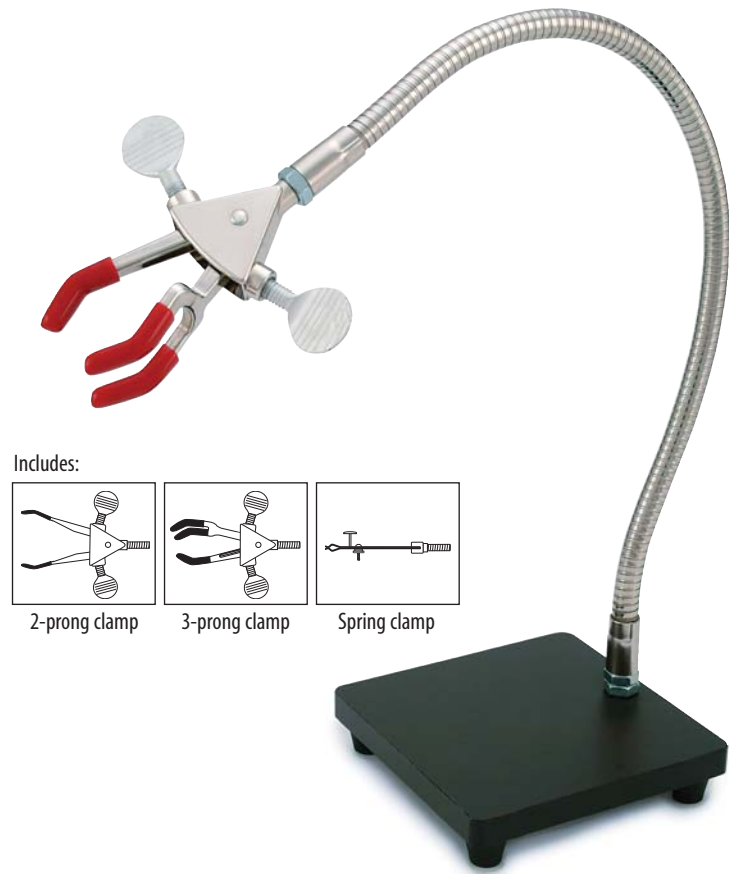


Ultra Flex Support System with Base Plate

- Stable base plate
- Unique flexible support arm
- All metal construction
- Includes 3 types of clamp heads

Ultra Flex Support System with Base Plate is designed with a stable base plate that fits easily on bench tops or in fume hoods. System includes the three interchangeable clamp heads that make it extremely versatile; 2-prong clamp head, 3-prong clamp head and spring clamp head. The system also includes the nickel plated flex arm in 12" (30.5cm) or 18" (45.7cm) length, steel base plate with chemical resistant black finish and utility wrench.

Specifications		
Arm Length	Ultra Flex 12	12" (30.5cm)
	Ultra Flex 18	18" (45.7cm)
Arm Diameter		0.5" (1.3cm)
Base Plate Dimensions		5" x 5" (12.7cm x 12.7cm)
Max. Grip Size	2-prong Clamp Head	2.6" (6.6cm)
	3-prong Clamp Head	2.3" (5.8cm)
	Spring Clamp Head	0.5" (1.3cm)
Ship Weight	Ultra Flex 12	5lbs (2.3kg)
	Ultra Flex 18	5lbs (2.3kg)

Ordering Information		
Description	Part Number	Price
Ultra Flex 12 with Base Plate	916310	\$179.00
Ultra Flex 18 with Base Plate	916311	215.00

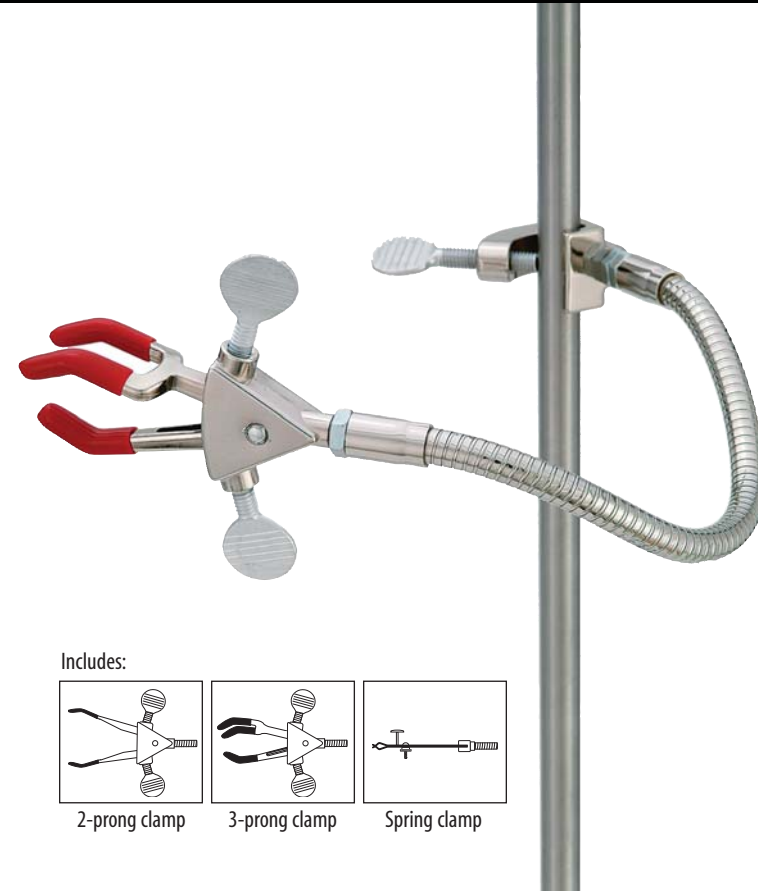


Ultra Flex Support System with Lab-frame Connector

- Mounts to lab-frames or support stands
- Ideal for use in fume hoods
- Unique flexible support arm
- Includes 3 types of clamp heads

Ultra Flex Support System with Lab-frame Connector securely attaches to support stands, lab-frames, or any support rod up to .75" (1.9cm) in diameter. Ultra Flex clamping system is extremely versatile. Unique, flexible arm can be placed in virtually any position or angle. System features 2-prong clamp head, 3-prong clamp head, spring clamp head, nickel plated flex arm in 12" (30.5cm) or 18" (45.7cm) length, cast alloy lab-frame connector and utility wrench.

Specifications		
Arm Length	Ultra Flex 12	12" (30.5cm)
	Ultra Flex 18	18" (45.7cm)
Arm Diameter		0.5" (1.3cm)
Max. Grip Size	2-prong Clamp Head	2.6" (6.6cm)
	3-prong Clamp Head	2.3" (5.8cm)
	Spring Clamp Head	0.5" (1.3cm)
	Lab-frame Connector	0.75" (1.9cm)
Ship Weight	Ultra Flex 12	2.1lbs (0.9kg)
	Ultra Flex 18	2.2lbs (1.0kg)

Ordering Information		
Description	Part Number	Price
Ultra Flex 12 with Lab-frame Connector	916314	\$139.00
Ultra Flex 18 with Lab-frame Connector	916315	165.00

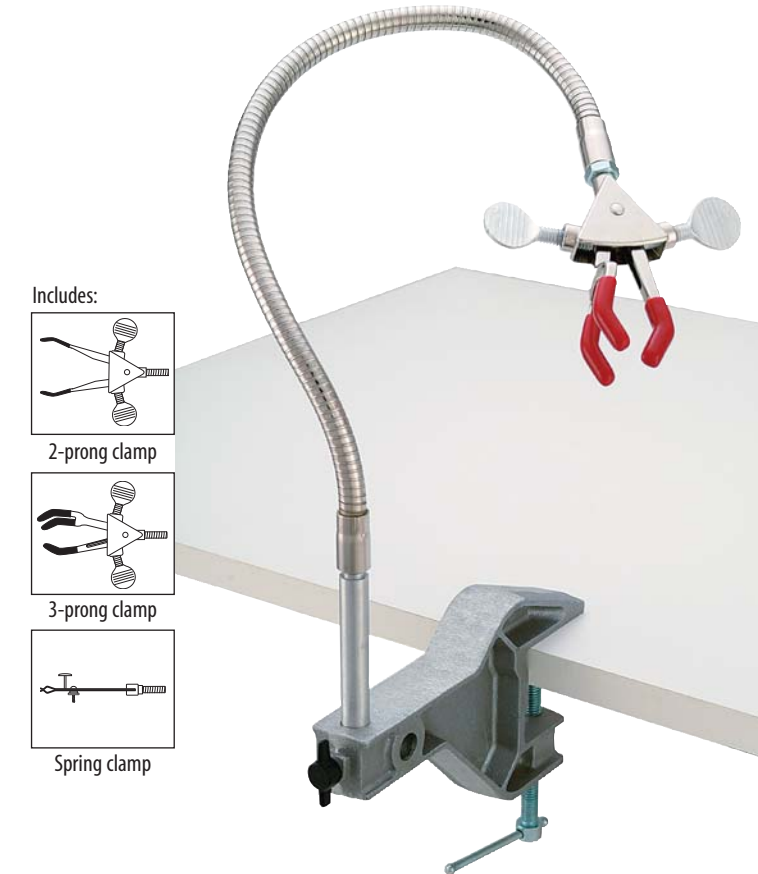

Ultra Flex Support System with Bench Clamp

- Easily attaches to most lab benches
- Frees up valuable bench space
- Unique flexible support arm
- Includes 3 types of clamp heads

Ultra Flex Support System with Bench Clamp easily attaches to the side of your lab bench or counter top, which helps utilize more space in your lab. Ultra Flex clamping system is extremely versatile with a unique, flexible arm that can be placed in virtually any position or angle. System includes 2-prong clamp head, 3-prong clamp head, spring clamp head, flex arm in 12" (30.5cm) or 18" (45.7cm) length, aluminum bench clamp and utility wrench.

Specifications		
Arm Length	Ultra Flex 12	12" (30.5cm)
	Ultra Flex 18	18" (45.7cm)
Arm Diameter		0.5" (1.3cm)
Max. Grip Size	2-prong Clamp Head	2.6" (6.6cm)
	3-prong Clamp Head	2.3" (5.8cm)
	Spring Clamp Head	0.5" (1.3cm)
	Bench Clamp	1.68" (4.3cm)
Ship Weight	Ultra Flex 12	2.4lbs (1.1kg)
	Ultra Flex 18	2.5lbs (1.2kg)

Ordering Information		
Description	Part Number	Price
Ultra Flex 12 with Bench Clamp	916312	\$159.00
Ultra Flex 18 with Bench Clamp	916313	190.00

All Ultra Flex Support Systems Include

3 types of interchangeable clamp heads designed to hold small apparatus such as PH meters, temperature probes, hydrometers, test tubes, vials, flasks and other small devices or containers.

2-prong Clamp Head



Constructed of cast alloy with dual adjust for greater flexibility. Ideal for holding straight-sided apparatus, such as, burets, thermometers, beakers and flask necks. Provided with vinyl and fiberglass sleeves. Fiberglass sleeves are recommended for higher temperature applications up to 140°C (284°F).

3-prong Clamp Head



More versatile and ideal for holding irregular shaped objects. Constructed of cast alloy with dual adjust for greater flexibility. Provided with vinyl and fiberglass sleeves. Fiberglass sleeves are recommended for higher temperature applications up to 140°C (284°F).

Spring Clamp Head



Features locking wing nut that applies tension to the nickel plated jaws. Designed to hold thermometers, glass tubes and temperature probes.

Replacement Parts

Description	Part Number	Price
2-prong Clamp Head	916321	\$18.00
3-prong Clamp Head	916322	24.00
Spring Clamp Head	916323	18.00
12" Flex Arm	916320	75.00
18" Flex Arm	916328	112.00
Base Plate	916324	79.00
Lab-frame Connector	916326	29.00
Bench Clamp	916325	49.00